




EZER™ Malaria Pf/Pan Antigen Rapid Test Device

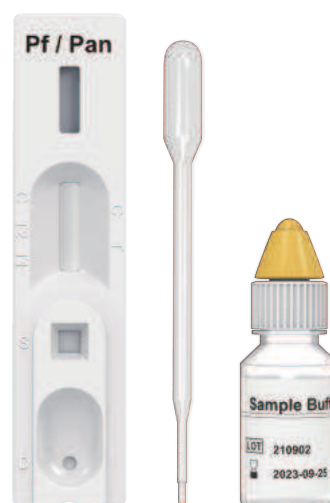

INTENDED USE

EZER™ Malaria Pf/Pan Antigen Rapid Test Device is a lateral flow chromatographic immunoassay for the qualitative detection of four kinds of circulating plasmodium (*P. falciparum*, *P. vivax*, *P. ovale*, and *P. malariae*) in human whole blood.

MALARIA



Hangzhou Genesis Biodetection and Biocontrol Co., Ltd.



+86-571 878 18163



www.genesis-ivd.com



310018



No.139, St.10th (East), Hangzhou Economic & Technological Development Zone.

PERFORMANCE CHARACTERISTICS

Result for Pf		Microscopy		
		+	-	Total
EZER™ Malaria Pf/Pan Antigen Rapid Test Device	+	47	0	47
	-	2	197	199
	Total	49	197	246


Relative Sensitivity:95.92% (95%CI*: 86.02% ~99.50%)
Relative Specificity:>99.99% (95%CI*: 98.14% ~100.00%)
Accuracy:99.19% (95%CI*: 97.09% ~99.90%)

Result for Pan		Microscopy		
		+	-	Total
EZER™ Malaria Pf/Pan Antigen Rapid Test Device	+	41	0	41
	-	1	197	198
	Total	42	197	239


Relative Sensitivity:97.62% (95%CI*: 87.43% ~99.94%)
Relative Specificity:>99.99% (95%CI*: 98.14% ~100.00%)
Accuracy:99.58% (95%CI*: 97.69% ~99.99%)

PROCEDURE

Whole Blood Specimen

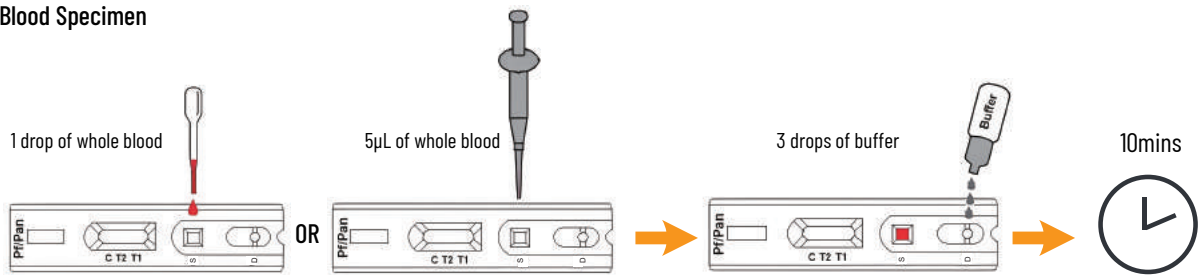

RESULTS

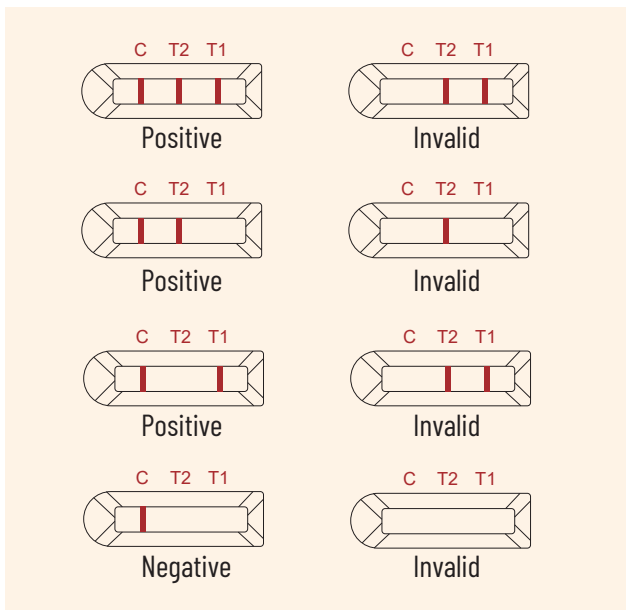

AVAILABILITY

Product	EZER™ Malaria Pf/Pan Antigen Rapid Test Device
Cat. No.	P231109
Test Device	20
Plastic Dropper	20
Sample Buffer	1
Package Insert	1