


EZER™ CTX-M β-lactamases Rapid Test



High Accuracy

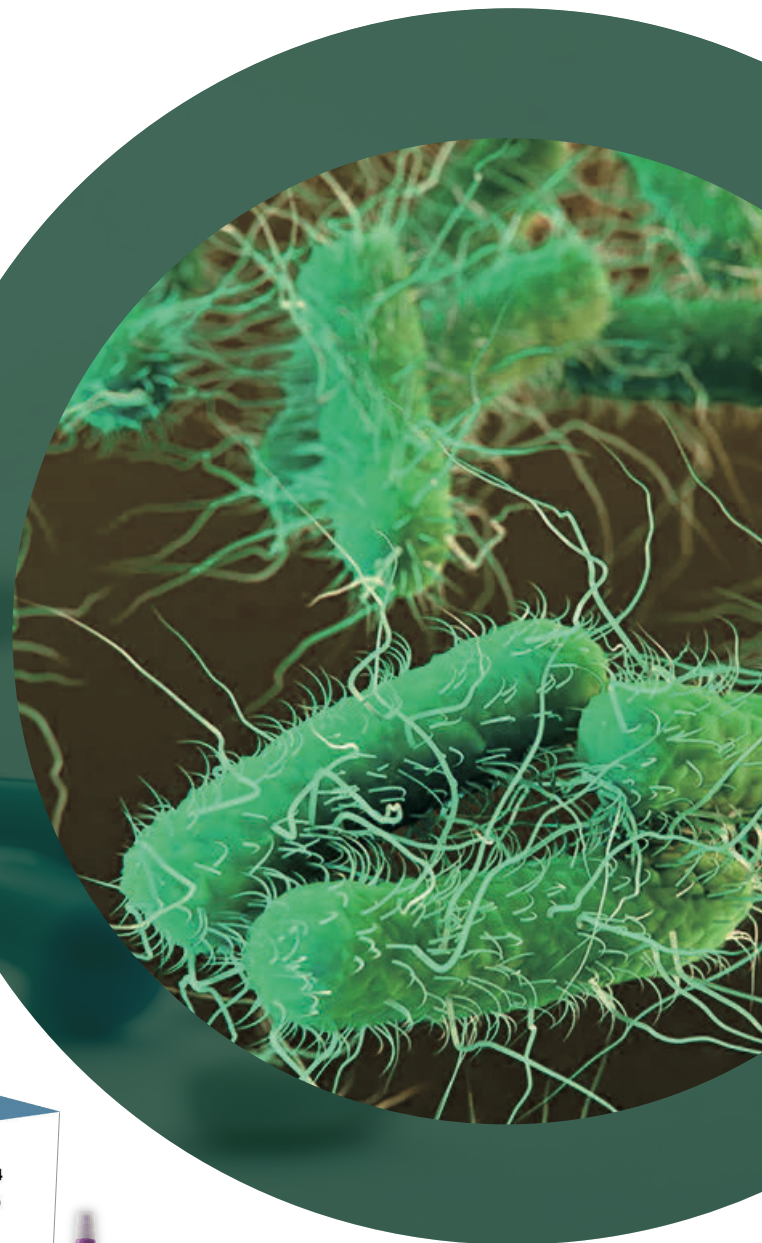

Fast Results in 15 Minutes



Easy to Use, No Equipment Required



Easy Visual Interpretation

INTENDED USE

The EZER™ CTX-M β-lactamases Rapid Test is a qualitative rapid immunoassay for the detection of CTX-M groups 1, 2, 8, 9 and 25 in a bacterial colony obtained from culture. It is an *in vitro* diagnostic assay, for professional use only, that aids to detect extended spectrum β-lactamases (ESBLs) produced by Enterobacterales.



+86-571 878 18163



310018



www.genesis-ivd.com



No.139, St.10th (East), Hangzhou Economic
& Technological Development Zone, Zhejiang, China


PERFORMANCE CHARACTERISTICS

Result		Other Commercial Test		
		+	-	Total
EZERT [™] CTX-M β -lactamases Rapid Test	+	135	0	135
	-	0	43	43
	Total	135	43	178

Relative Sensitivity: 100.00%
 Relative Specificity: 100.00%
 Accuracy: 100.00%


PROCEDURE

1 Bacterial Colonies

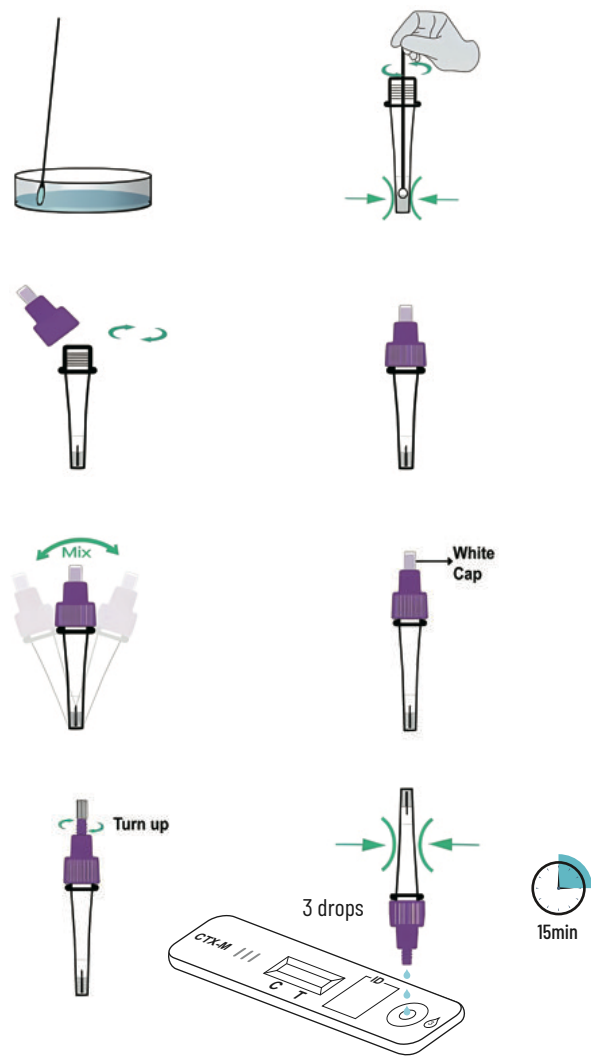

2 For positive blood culture samples, please use the *EZER[™]* BloodPrime for pretreatment.

RESULTS

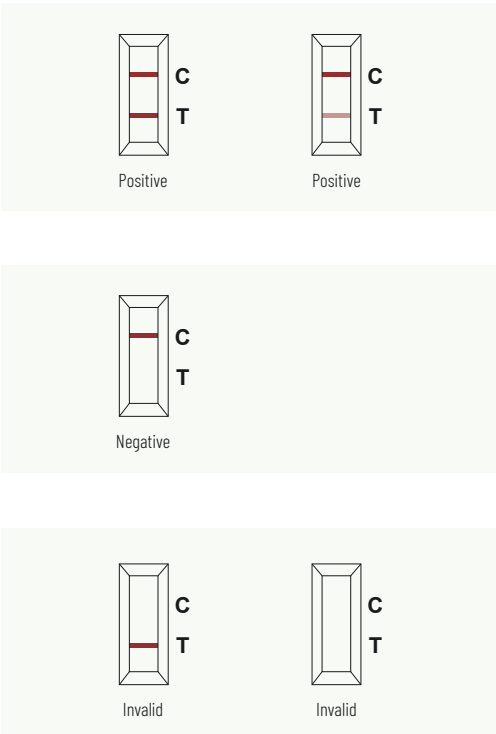

AVAILABILITY

Product	<i>EZER[™]</i> CTX-M β -lactamases Rapid Test
Cat. No.	P261104
Test Device	20
Extraction Tube with Buffer	20
Package Insert	1