


EZER™ OXA-23 Carbapenemase Rapid Test



High Accuracy Compared with PCR

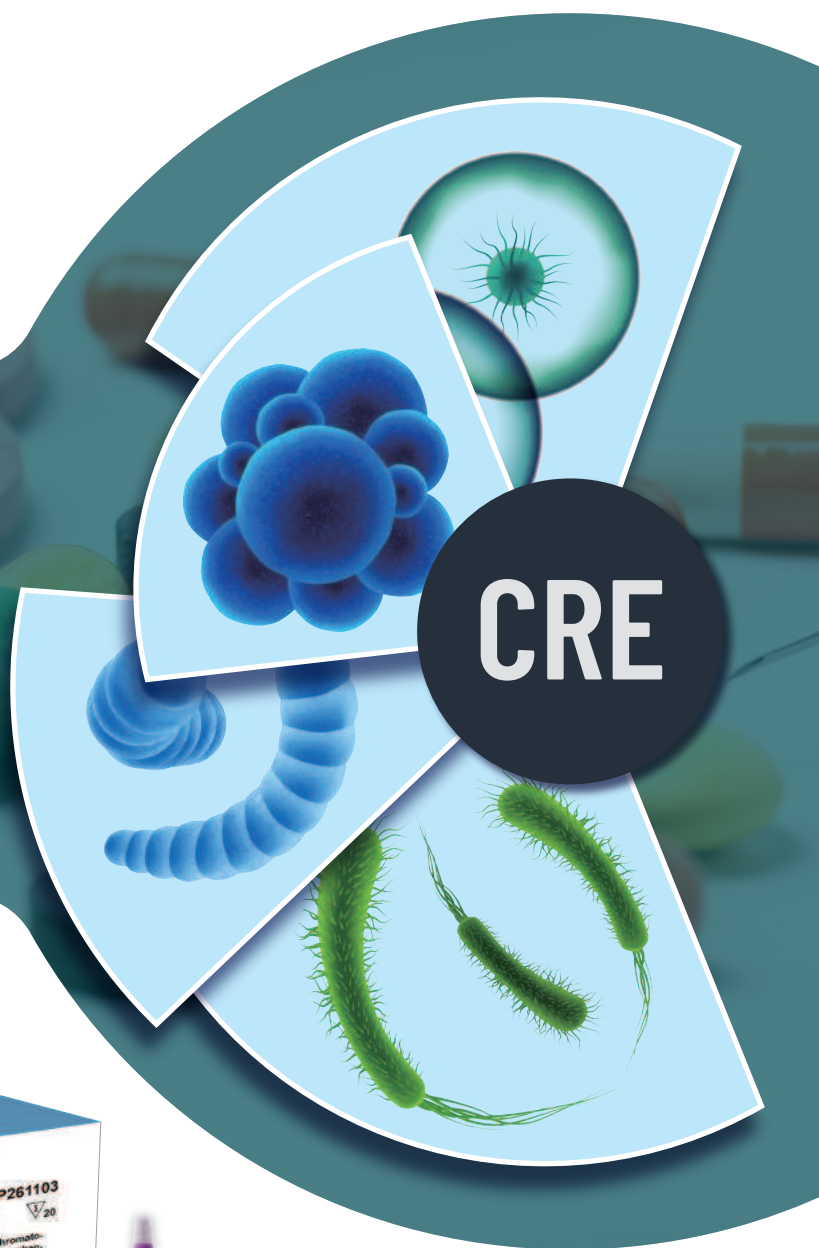

Fast Results in 15 Minutes



Easy to Use, No Equipment Required



Easy Visual Interpretation

INTENDED USE

The EZER™ OXA-23 Carbapenemase Rapid Test is an immunochromatographic assay (ICA) for directly qualitative detection of OXA-23 carbapenemases in bacterial culture. It's an in vitro diagnostic test for professional use only.



+86-571 878 18163



310018



www.genesis-ivd.com




No.139, St.10th (East), Hangzhou Economic
& Technological Development Zone, Zhejiang, China

PERFORMANCE CHARACTERISTICS


Result		Molecular Method		
		+	-	Total
EZER™ OXA-23 Carbapenemase Rapid Test	+	235	0	235
	-	1	38	38
	Total	235	38	273

Relative Sensitivity: 100.00%
 Relative Specificity: 100.00%
 Accuracy: 100.00%

PROCEDURE

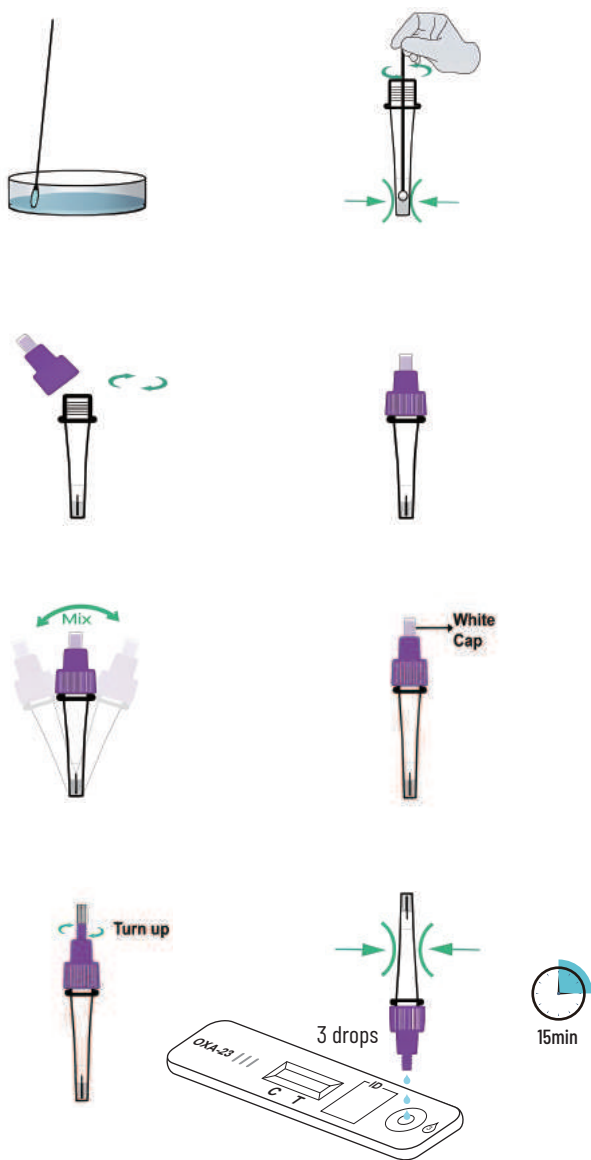

RESULTS

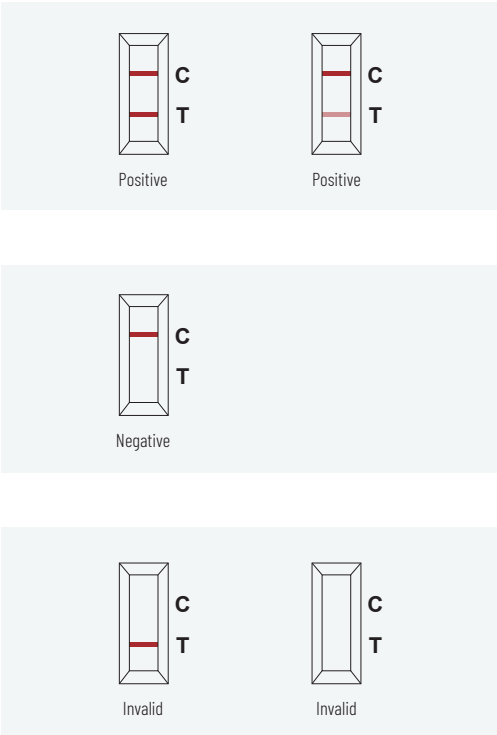

AVAILABILITY

Product	EZER™ OXA-23 Carbapenemase Rapid Test
Cat. No.	P261103
Test Device	20
Extraction Tube with Buffer	20
Package Insert	1